

FMT

Medical Sensors, Cables and Accessories



PRESSURE INFUSION CUFF



Applications :

- Surgery • Intensive Care • Radiology, Special Procedures • Cardiac Catheterization Laboratory • Emergency/Trauma Department

FMT Reusable Pressure Infusor is used for the rapid liquid infusion, in invasive pressure monitoring and autotransfusion procedures.

FMT Reusable Pressure Infusor is made of thick, waterproof polyester fabric provides most durable quality, delivers the maximum pressure and assures cleanness. The rear side of the inflation bag is supported with internal solid plates to prevent the excessive deformation of the pressure infusor during inflation.

Polyester band I.V. pole loop hanger and durable polyamide I.V. fluid bag hook are designed to provide long service life. Fluid levels are easily viewed through the transparent net sleeve. Accurate pressure readings are displayed on precision aneroid gauge. FMT pressure infusor comes with a palm fitted efficient hand bulb with thumbwheel air valve which allows single hand operation.

FMT pressures infusor are adaptable with most I.V. bags and available in 500ml, 1000ml, and 3000ml sizes.

Features:

- Durable, polyester band I.V. Pole loop hanger
- Accurate pressure gauge is easy to read and durable, ± 2 mmHg at 300mmHg.
- I.V. Fluid Bag Hook is durable and designed to fit all Standard bags snugly.
- Durable Construction - will not rip or tear.
- Clear sleeve around bag provides a visible window to see what fluids are being pressurized.
- Palm fitted, efficient hand bulb makes inflation easy.
- TPU Bladder
- Washable up to 30°C.
- Latex free.



- IB 3000 3000 ml. replacement bladder.
- IB 1000 1000 ml. replacement bladder.
- IB 500 500 ml. replacement bladder.
- PI MANO Replacement Manometer.
- PI BULB Replacement Hand Bulb.

PIHP30 Hand Pump for Pressure Infusors
(For use with 3000 ml Infusors Only)



FMT-PI/05 Reusable Pressure Infusor, 500 ml

Pocket type, with hand bulb and manometer.



FMT-PI/10 Reusable Pressure Infusor, 1000 ml

Wrap-Around type with velcro, with hand bulb and manometer.



FMT-PI/10S Reusable Pressure Infusor, 1000 ml

Pocket type, with hand bulb and manometer.



FMT-PI/30 Reusable Pressure Infusor, 3000 ml

Pocket with zip type, with hand bulb and manometer.



Applications:

- Surgery
- Intensive Care
- Radiology, Special Procedures
- Cardiac Catheterization Laboratory
- Emergency/Trauma Department

FMT Disposable Pressure Infusor is used for the rapid liquid infusion, in invasive pressure monitoring and autotransfusion procedures. The FMT Disposable Pressure Infusor is designed for single patient use potentially preventing and controlling cross-contamination and hospital acquired infections.

FMT Disposable Pressure Infusor is made of polyurethane laminated fabric provides most durable quality, delivers the maximum pressure and assures cleanness. Fluid levels are easily viewed through the mesh sleeve. The safety valve eliminates the possibility of over inflation and releases air at 330 mm/Hg. FMT pressure infusor comes with a palm fitted efficient hand bulb with 3-way stopcock air valve.

FMT pressures infusors are adaptable with most IV bags and available in 500ml and 1000ml sizes.

FMT pressure infusors are Latex Free.

Features:


- Designed for single patient use.
- Safety valve prevents over-inflation (330 mmHg pressure relief).
- Color-coded gauge makes monitoring pressure simple (0-300 mmHg).
- 3-way stopcock ensures precise control of pressure.
- Durable, metal sleeved IV Pole Loop Hanger.
- IV fluid bag strip is durable and designed to fit all Standard bags snugly.
- Mesh sleeve around bag provides a visible window to see what fluids are being pressurized.
- Palm fitted, efficient hand bulb makes inflation easy.
- Latex free.



FMT-DPI/05
FMT-DPI/10

Disposable Pressure Infusor, 500 ml
Disposable Pressure Infusor, 1000 ml

FMT is registered trademark of **Metko Ltd.**

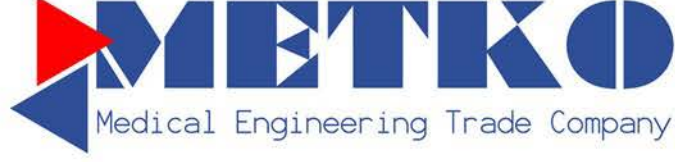

Headquarter & Factory

İvedik O.S.B. Ağaç İşleri Sanayi Sit.
1354. Cad. 1358.Sok. No:9
06378 Yenimahalle

ANKARA / TÜRKİYE

Tel (pbx) : +90 312 387 12 46
Fax : +90 312 387 12 51

Branch

Atatürk Mah. Sedef Cad. 38.Ada
Ata Plaza 3-5 Blokları Kat:2 No:57
34758 Ataşehir

İSTANBUL / TÜRKİYE

Tel (pbx) : +90 216 455 91 21
Fax : +90 216 455 48 24

www.metkold.com

Other products, logos and company names mentioned herein may be trademarks of their respective owners and are not to be taken as an endorsement of affiliation with FMT®. The information in this catalogue has been carefully checked and it is believed to be accurate. Images in this catalogue are representative of items offered. Actual product may have slight visual differences. In the interest of continued product development, Metko reserves the right to make changes and improvements to this catalogue and the product it described any time, without notice or obligation.